



DR. BEN LARRABEE, D.D.S., M.S.


Dr. Larrabee has been practicing orthodontics for over 10 years now, with all of his years being in the competitive markets of Chicago and Phoenix. As a graduate from the University Michigan Dental School, Marquette Orthodontics, and friend and mentorship of Dr. David Sarver - progressive orthodontic treatment has always been a passion. However, always being the new practice in a competitive area has also created another passion - creating the presence of his care and practice in the community he serves. With the wonderful teachings from Dr. Terry Sellke to LeAnne Peniche, Dr. Larrabee has helped create a progressive presence for his area - and hopes to help you to create yours.

AREAS OF INTEREST: SELF LIGATION | DUAL ACTIVATION | EMPOWER | PRACTICE MANAGEMENT

Empower[®] Self Ligating Braces

Empower is the industry's most complete self ligating bracket system, bringing you choice like you've never had before. This unprecedented array of possibilities inspires you to treat patients your way. Empower puts control back in your hands.

- Unique design options allow you to adapt your treatment philosophies to meet your patients' needs
- Low profile, patient-friendly hooks (optional)
- Empower Clear - Reliable bond strength of proven Quad Matte™ base exclusive to American Orthodontics
- Empower SL - Maximum Retention pad with photo-etched pockets beneath 80 gauge mesh to increase pad surface area and ensure safe, secure bonds

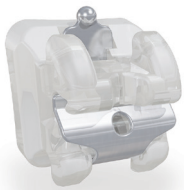

EMPOWER CLEAR

Clear bracket combines performance with patient-pleasing aesthetics

- Durable, aesthetically plated clip
- Individual color coding and Visual Positioning Aids
- Chamfered slot entrances to reduce wire binding
- Reliable bond strength of proven Quad Matte™ base exclusive to AO



EMPOWER METAL

Metal brackets offer fully interactive, fully passive, or Dr. Larrabee's choice, combination Dual Activation™ options

- Precision heat treated chromium cobalt clip
- Extended tie wings provide easy tie-ability and allow differential tying
- Center line allows accurate bonding alignment with long axis of tooth



ifit MOLAR TUBES

The ultimate combination of functionality, convenience, and patient comfort.


- Near horizontal hook orientation
- Enhanced bond strength and fit
- Ultra low profile

Empower SL

DUAL Activation[™] SYSTEM

The Dual Activation System is a single, innovative concept with multiple applications - all designed to enhance your treatment philosophy. Here's what you'll experience with Dual Activation:

- The best of both worlds - interactive and passive brackets in one unified system
- Excellent anterior control
- Improved posterior freedom play
- Time saving benefits and ease of self ligation

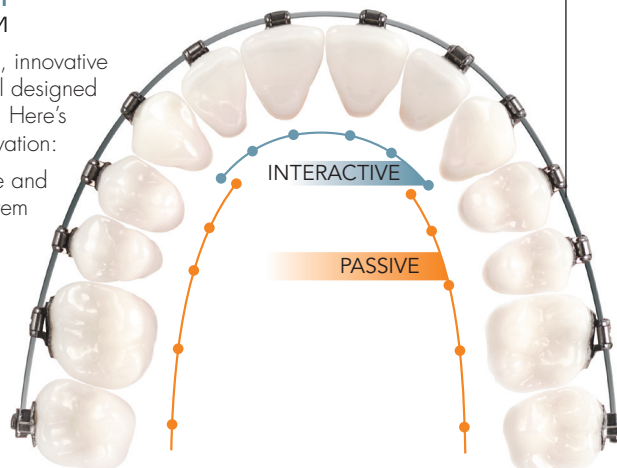

PERIPHERALS

LUNO SUSHI PLIERS 001-017403

Ceramic Bracket Debonding

- American forged steel
- Sleek two-tone finish
- Smooth, flush joint



EMPOWER INSTRUMENT 678-901C

Opening Instrument

- Single prong tip for opening using gingival tail of clip
- Double prong tip for opening using facial hole or for pressing wire into bottom of slot during clip closing



DEBONDING BURS

Carbide Debonding Burs

- Controlled debonding without nicking, scratching, or abrading the enamel



TANZO WIRE

Copper Nickel Titanium Wire

- Polished surface finish
- Preinstalled stops - optional
- High resiliency to deformation
- Quality ensured



TANZO

ELASTOMERICS

Wildlife Elastics

- Top quality latex tubing
- Precision cutting ensures consistent pulling forces
- Latex, Latex Free, and Neon options



DR. LARRABEE'S RECOMMENDED EMPOWER DUAL ACTIVATION PRESCRIPTION



Empower
Self Ligating

.022

Maxillary	Torq	Ang	Rot	M-D	R/L	INTERACTIVE	PASSIVE
Central	+12	+5		Empower SL	R L	485-1112 485-2112	
Lateral	+8	+9		Empower SL	R L	485-128 485-228	
Cuspid	0	+8		Empower SL With Hook	R L	485-130B 485-230B	
1st Bicuspid	-7	0		Empower SL - Standard Pad With Hook	R L		585-1407B 585-2407B
2nd Bicuspid	-7	0		Empower SL - Standard Pad With Hook	R L		585-1407B 585-2407B
1st Molar	-18		10	ifit Direct Bond Tube Non-Convertible, With Hook	4.5 mm R L		1957-22-000047HM 1958-22-000048HM
2nd Molar	-25		5	ifit Direct Bond Tube Non-Convertible, With Hook	3.2 mm R L		1951-22-016220HM 1952-22-016221HM
Mandibular							
Anteriors	-1	0		Empower SL	R L	485-4101 485-3101	
Cuspid	0	+3		Empower SL With Hook	R L	485-430B 485-330B	
1st Bicuspid	-12	+2		Empower SL - True Offset Pad With Hook	R L		585-4412B 585-3412B
2nd Bicuspid	-17	+2		Empower SL - True Offset Pad With Hook	R L		585-4517B 585-3517B
1st Molar	-28		0	ifit Direct Bond Tube Non-Convertible, With Hook	4.5 mm R L		1955-22-000035HM 1956-22-000036HM
2nd Molar	-10		5	ifit Direct Bond Tube Non-Convertible, With Hook	3.2 mm R L		1953-22-016228HM 1954-22-016229HM



Empower
Clear Self Ligating

.022

Maxillary	Torq	Ang		R/L	INTERACTIVE	PASSIVE	
Central	+17	+4		Empower Clear	R L	1485-22-1117 1485-22-2117	
Lateral	+10	+8		Empower Clear	R L	1485-22-1210 1485-22-2210	
Cuspid	0	+8		Empower Clear With Hook	R L	1485-22-130HD 1485-22-230HD	
1st Bicuspid	-7	0		Empower Clear With Hook	R L		1585-22-1407HD 1585-22-2407HD
2nd Bicuspid	-7	0		Empower Clear With Hook	R L		1585-22-1407HD 1585-22-2407HD



EMPOWER MOLAR TUBES

Reconvertible molar tubes provide versatility without compromise

- Multiple design choices - single or double, direct bond or weldable
- Passive clip design for low friction sliding mechanics and easy wire insertion
- Full slot-width clip delivers exceptional rotation control

